**SPDC Data Collection
Instruction Manual**

**SPDC Data Collection
Manual de instrucciones**

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1. INTRODUCTION

This software automatically detects data from the instrument and enables the user to record the received data into a Text, Word, Excel, Access or PDF file. PDF file is for Starter/Aquasearcher series only.

2. INSTALLATION

Execute the “setup.exe” file to install the software.


3. START-UP

After installation has been completed, start the software by clicking on either of the following:

- The “SPDC DataCollection” icon on desktop.
- The “SPDC DataCollection” icon in programs menu.
- The “SPDC DataCollection” file in installation folder.

4. MAIN WINDOW


When the software has been launched, the main window will open.

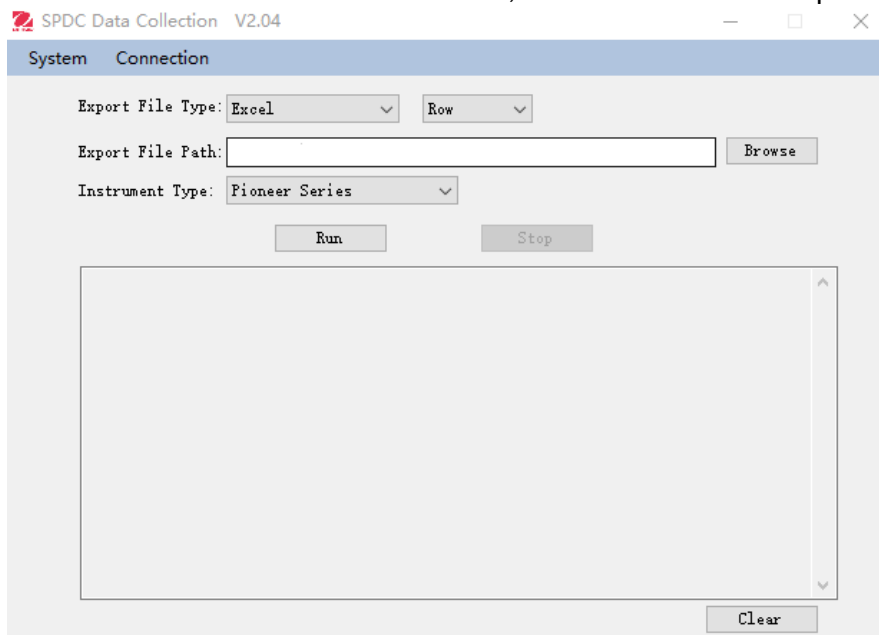5. SERIAL PORT CONFIGURATION

Click the “Config” link in “Connection” menu, the <Connection Configuration> window will appear.

Click the “RS232”. The default configuration is 9600-8-N-1-N.

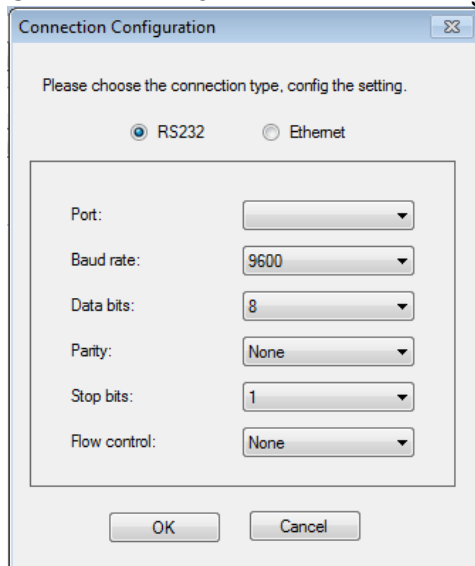The screenshot shows the 'Connection Configuration' dialog box with the 'RS232' radio button selected. The settings are as follows:

Parameter	Value
Port:	[Dropdown menu]
Baud rate:	9600
Data bits:	8
Parity:	None
Stop bits:	1
Flow control:	None

Buttons: OK, Cancel

The following can be set:


- The serial port number to which the instrument is connected.
- The baud rate.
- The number of data bits.
- The parity.
- The number of stop bits.
- The flow control.

Warning: The baud rate, the number of bits, the parity and the number of stop bits selected in the software must correspond exactly to those of the instrument settings. Otherwise, the data will not be received or will be misinterpreted.

6. ETHERNET CONFIGURATION

Click the “Config” link in “Connection” menu, the <Connection Configuration> window will appear.

Click the “Ethernet”. Then the IP address and Port can be set.

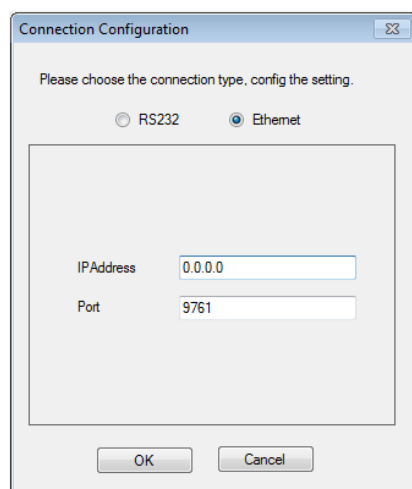The screenshot shows the 'Connection Configuration' dialog box with the 'Ethernet' radio button selected. The settings are as follows:

Parameter	Value
IPAddress	0.0.0.0
Port	9761


Buttons: OK, Cancel

Warning: The IP address and Port filling in the software must correspond exactly to those of the instrument settings. Otherwise, the data will not be received or will be misinterpreted.

7. EXPORT FILE TYPE

Click the dropdown list, and choose the type of file that you want to export the received data to from the instrument.

- Text File
- Word File
- Excel File
- Access File
- PDF File (when selecting "Starter/Aquasearcher Series" as Instrument Type)

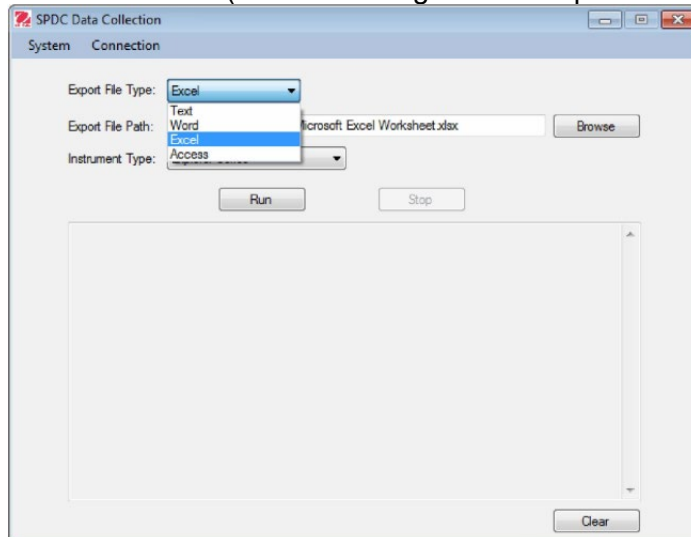
8. EXPORT FILE PATH

According to the Export File Type, you may choose a file to which you want to export the received data.

- Input file path manually or choose file by clicking "Browse" button.
- The export file can be an existing file, or you may create a new file.

9. EXPORT INSTRUMENT TYPE

Click the dropdown list and choose the type of instrument that you want to link.

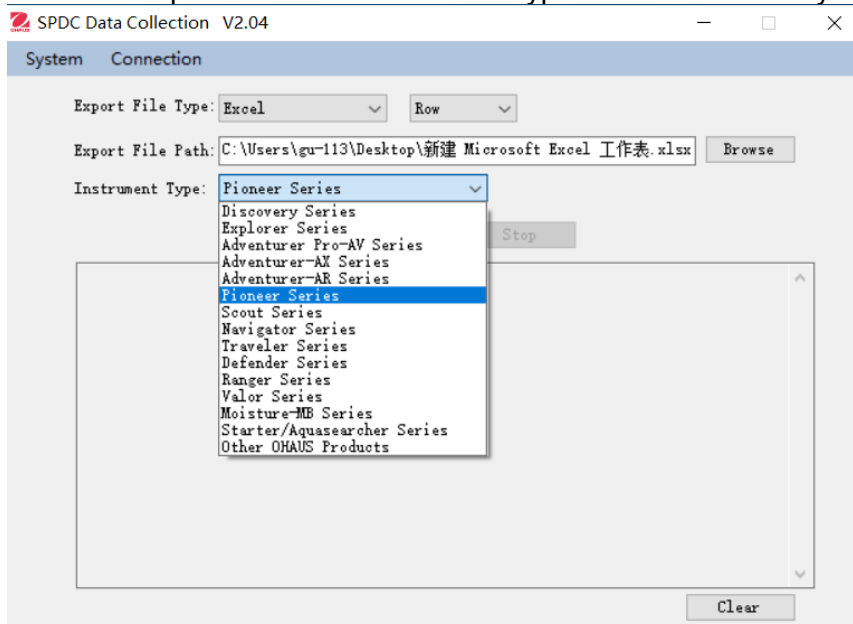Data Collection

1. Prerequisite:

- The serial port settings must correspond exactly to instrument settings.
- The chosen export file must exist and be writable.

2. Start:

Click the “Run” button.

The software will now detect the data received from the balance.

When there are data detected from instrument, the software will export the received data to both file and screen automatically.

Note: For Word and Excel files, the file will be opened and the received data will display synchronously.

3. Field input:

For Excel files, the input field could be customized by point the cursor to certain grid.

4. Stop:

Click the “Stop” button.


The software will stop detecting data received from the instrument.

Note: For Word and Excel files, if the file is closed the detection will also stop.

5. Clear data:


Clear data in the screen by clicking the “Clear” button.

Examples: Data transfer to Word file.

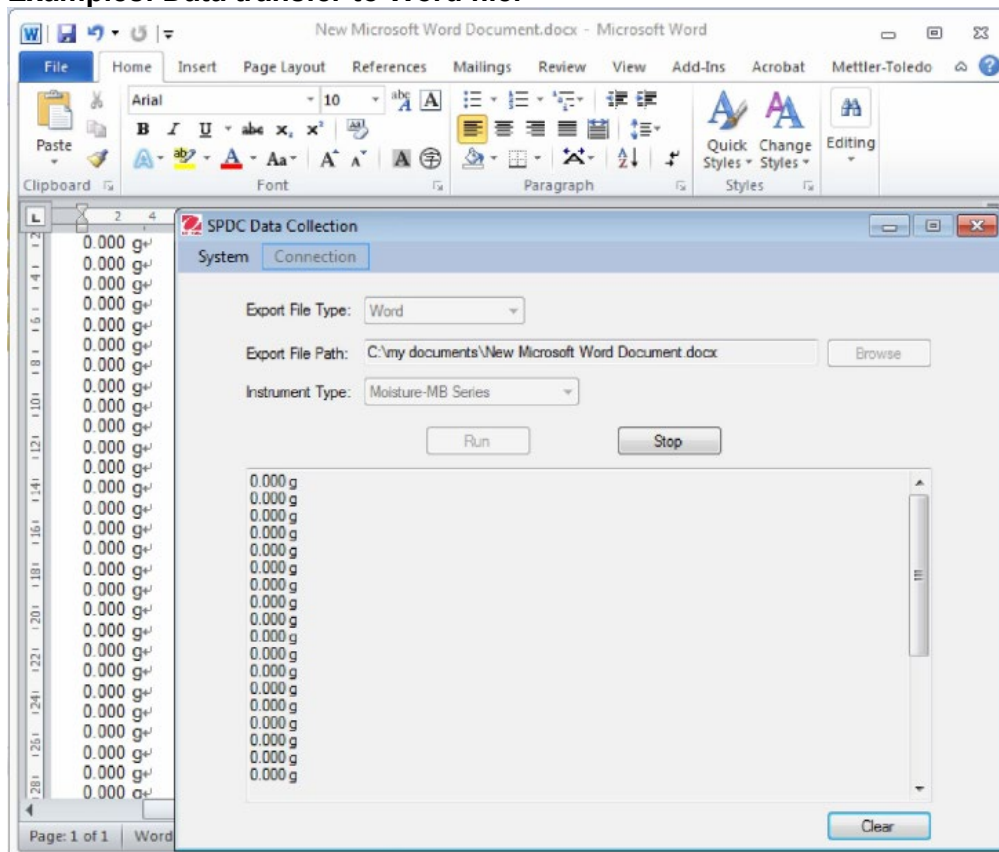
Examples: Data transfer to Excel file.

- In vertical format (when selecting “Row”)

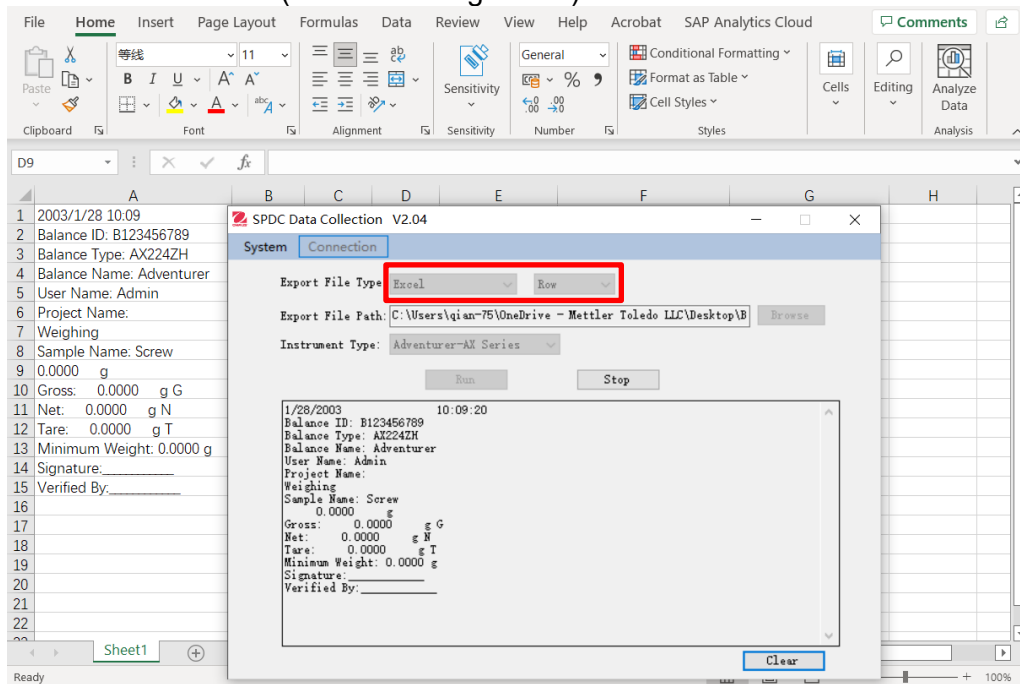
- In horizontal format (when selecting “Column”)

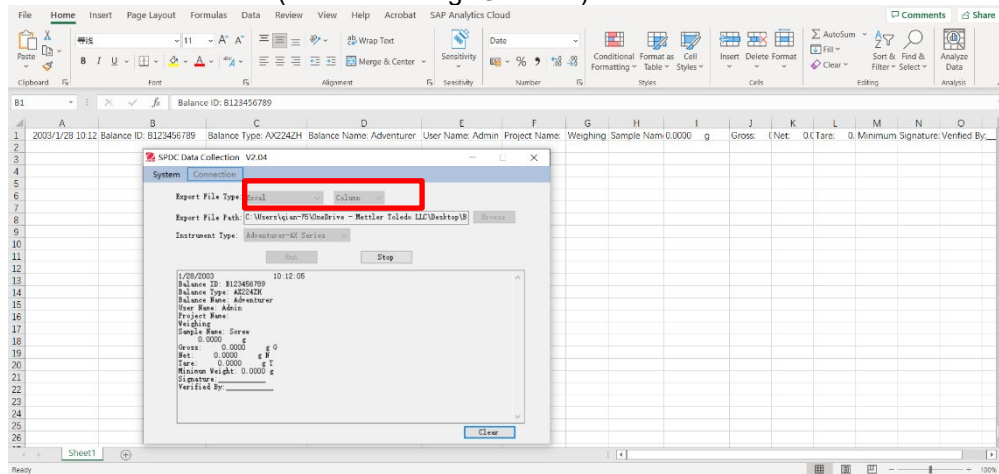Example: Data transfer to PDF file (for Starter/Aquasearcher Series only)

Please select “Starter/Aquasearcher Series” when connecting AB33EC, AB33M1, AB33PH and AB41PH meters to PC.

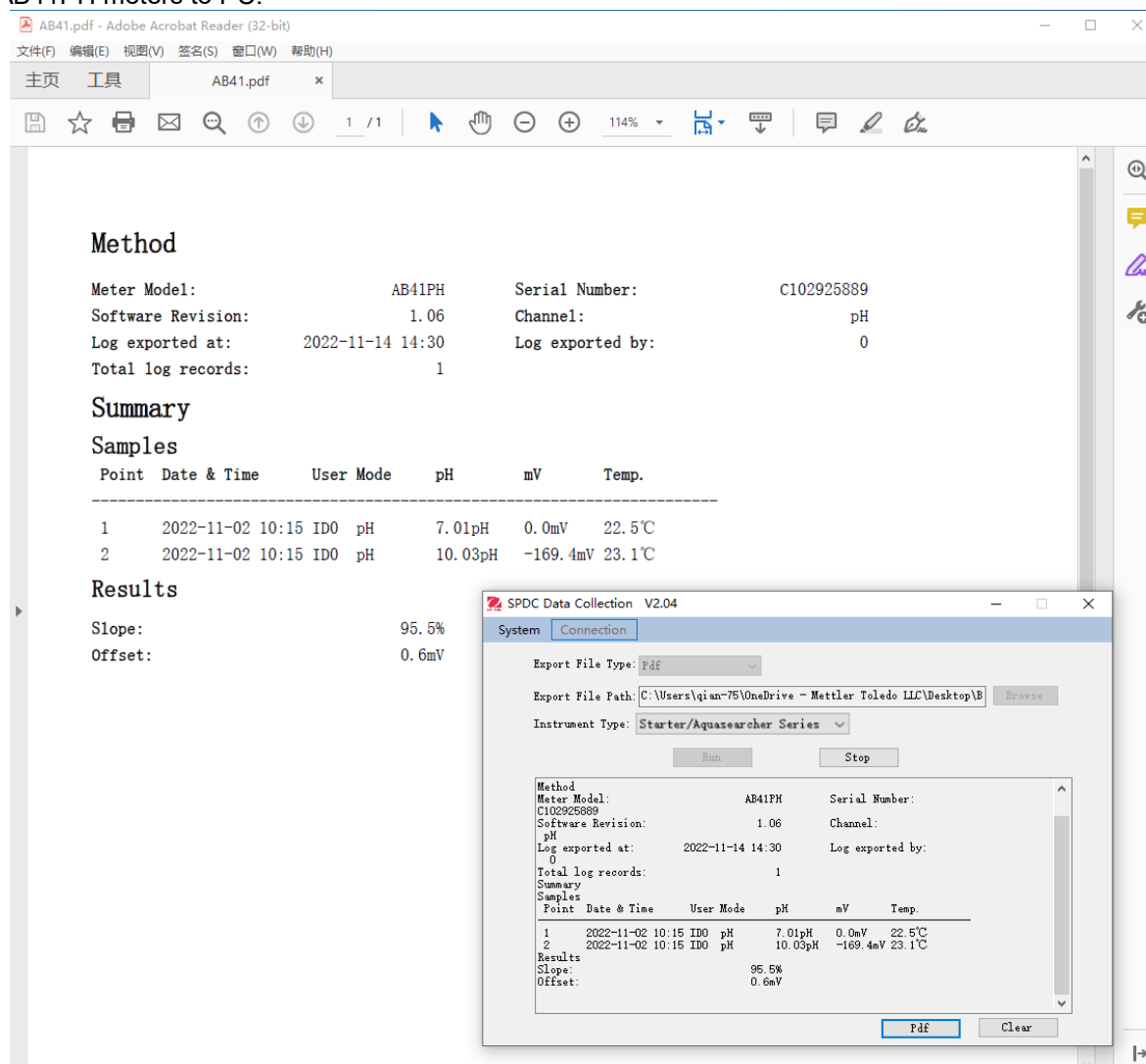**10. EXIT**

Exit the software by clicking:

- _ The “Exit” link in “System” menu.
- _ The “X” button in software main window.

11. UNINSTALL

Uninstall the software by clicking:

- _ The “Uninstall” icon in programs menu.
- _ The “Uninstall.exe” file in installation folder.

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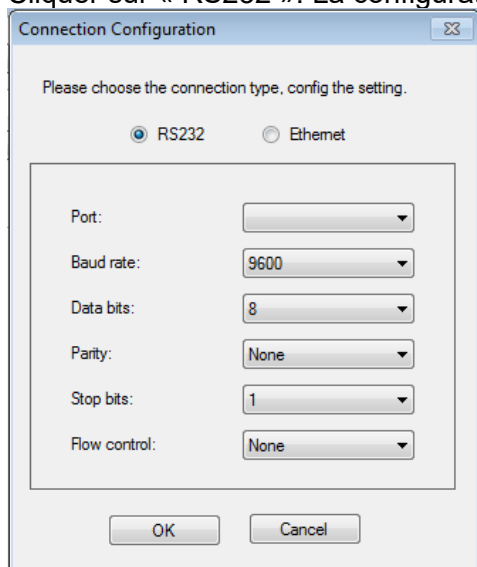
The screenshot shows a 'Connection Configuration' dialog box. At the top, it says 'Please choose the connection type, config the setting.' Below this, there are two radio buttons: 'RS232' (which is selected) and 'Ethernet'. Underneath, there is a list of settings, each with a dropdown menu: 'Port:', 'Baud rate:' (set to 9600), 'Data bits:' (set to 8), 'Parity:' (set to None), 'Stop bits:' (set to 1), and 'Flow control:' (set to None). At the bottom of the dialog are 'OK' and 'Cancel' buttons.

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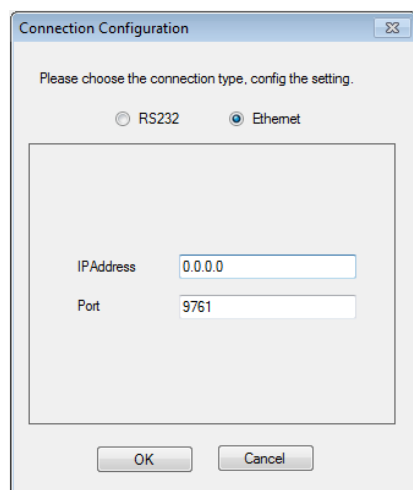
The screenshot shows the 'Connection Configuration' dialog box with 'Ethernet' selected. The 'RS232' radio button is unselected, and the 'Ethernet' radio button is selected. The settings area now contains two text input fields: 'IPAddress' with the value '0.0.0.0' and 'Port' with the value '9761'. The 'OK' and 'Cancel' buttons are at the bottom.

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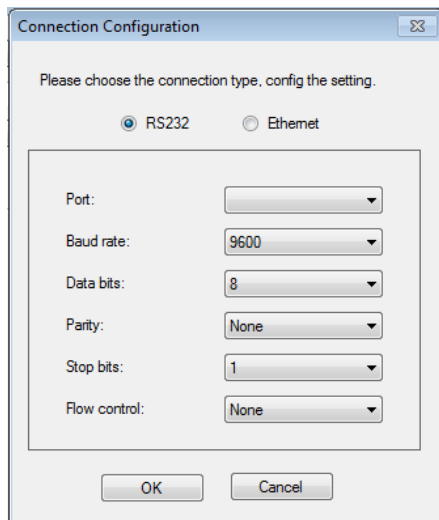
The dialog box is titled "Connection Configuration" and contains the instruction "Please choose the connection type, config the setting." Below this, there are two radio buttons: "RS232" (selected) and "Ethernet". Underneath, there is a group box containing several settings, each with a label and a dropdown menu: "Port:" (empty), "Baud rate:" (9600), "Data bits:" (8), "Parity:" (None), "Stop bits:" (1), and "Flow control:" (None). At the bottom of the dialog are "OK" and "Cancel" buttons.

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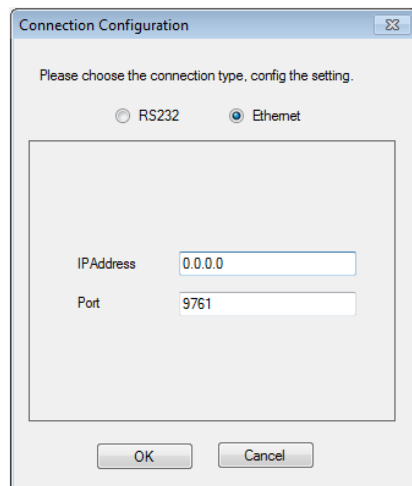
The dialog box is titled "Connection Configuration" and contains the instruction "Please choose the connection type, config the setting." Below this, there are two radio buttons: "RS232" and "Ethernet" (selected). Underneath, there is a group box containing two settings, each with a label and a text input field: "IPAddress" (0.0.0.0) and "Port" (9761). At the bottom of the dialog are "OK" and "Cancel" buttons.

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P/N 30240594 J © 2022 OHAUS Corporation, all rights reserved / todos los derechos reservados / tous droits réservés / Alle Rechte vorbehalten / tutti i diritti riservati
Ohaus Corporation